




Wheel Lock Extension

If you have any questions call Sunrise Medical Technical support at 800-333-4000 or visit www.SunriseMedical.com

A. WHEEL LOCK EXTENSION REPLACEMENT KIT:

Remove existing nut and discard (A). Replace extension handle with style ordered (B or C) and install with new nut included in the kit. Tighten until mechanism does not move freely then back off 1/4 turn.

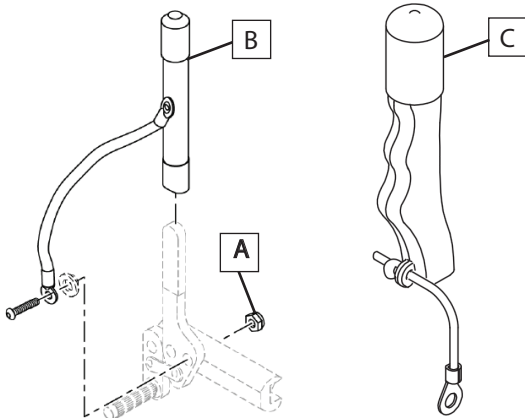

Sunrise Medical 2842 Business Park Ave
Fresno, CA 93727 USA
In Canada (800) 263-3390



© Sunrise Medical (US) LLC
06.17 245687 Rev. B